hairegen

User Manual

Ministry of Health registration 30870003

Distributed by Dimedico Medical Technologies Ltd.

91 Sokolov St. Ramat Hasharon

Tel: 03-5255722 | Fax: 03-5255822

info@dimedico.co.il | www.dimedico.co.il

Manufactured by Pilogics LLP

30 Ha-Kharoshet St. Or Yehuda

Tel: 03-3731382 | Fax: 03-6911686

INDEX

- Warnings
- About Hairegen
- Using Hairegen
- Maintenance

WARNINGS

- Please read the instructions before use.
- Keep the device away from children.
- The device is not to be used by children. Persons with reduced physical, sensory or mental capabilities may use Hairegen only under supervision.
- "Class 1" laser product. Do not shine the red laser light into the eyes*.
- Use on the scalp only - do not use on other body parts.
- Do not remove the cartridge while the applicator is turned on.
- Keep the device away from water or wet areas.
- Do not use on wet hair.
- Store and operate the device at normal room temperature.
- Avoid using the device if it was dropped, physically damaged or wet. In case of malfunction please contact the seller, do not try to fix the problem by yourself.
- Do not use the device if the cartridge discs are damaged.

Do not use Hairegen in case of:

- History of contact allergies.
- Pregnancy.
- Scars, moles, wounds or other lesions on the treated scalp areas.
- Scalp disorders such as eczema, seborrheic dermatitis, psoriasis, actinic (solar) keratosis, acne, rosacea, infections, other inflammatory conditions or other chronic conditions.
- Treatment now or within the past 6 months with retinoids.

Discontinue using Hairegen if experiencing any significant discomfort during or after usage. Excessive use beyond recommendations may cause discomfort or irritation.

* In some models

ABOUT HAIREGEN

Hairegen is a personal use scalp stimulation device that helps improve hair thickness and density. When used according to instructions, Hairegen provides the scalp with multiple forms of stimulations: optical *, physical, electrical and biochemical. The device is suitable for men and women.

HAIREGEN KIT CONTAINS

1. Applicator

The battery operated applicator powers the cartridge.

2. Replaceable cartridge

The cartridge is the component that comes in contact with the scalp and delivers the stimulations. The cartridge must be connected to the applicator. The cartridge is for personal use only. The cartridge has limited usage time and can be replaced by the user.


3. Base

The base is a docking and charging station for the applicator. The base also controls the applicator's stimulation profile and enables communication with a computer.

4. USB cable and power adaptor

A standard micro-USB cable provides power to the base and data connectivity to an optional PC or Mac data presentation software (Personal Trainer). The USB cable connects to a standard USB power adapter or a powered USB port.

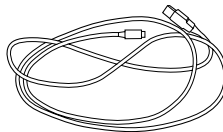
PREPARING FOR USE

1. Connecting the USB cable to the base


Connect the USB cable to the port at the bottom of the base. Connect the other end of the USB cable to the provided power supply, or a powered USB port.

Place the base on a stable flat surface, with the cable facing down and the applicator's cradle facing up.

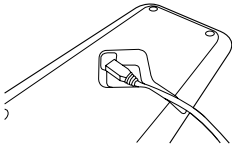2. Attaching the cartridge to the applicator

Remove the cartridge from its packaging. Ensure that the disks roll freely by spinning them a couple of times. Verify that the disks are not bent. Firmly push the cartridge into the applicator until it clicks into place.

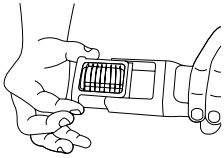
3. Charging the applicator

The base should be connected to a USB power source via the cable. The cartridge should be firmly connected to the applicator. The applicator should be fully inserted into the base. If unsure, remove and replace.

When cradled in the base, the applicator's light turns on, and the configuration buttons on the base turn on a few seconds later. If the battery needs charging, the applicator's light will be orange while charging.


When the battery is fully charged, the light turns green. It is recommended to keep the applicator in the base and have the base powered at all times.

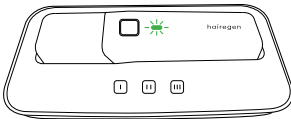4. Selecting a stimulation profile

When the base is powered and the applicator is properly cradled, the 3 touch buttons I - II - III are lit.

Each button corresponds to a specific stimulation profile. The button represents the chosen profile of the applicator. Profile selection is enabled only while the base is powered and the applicator is cradled. Touch a dimly lit button to select its corresponding stimulation profile. When selected, the active profile button will light up brighter than the others.

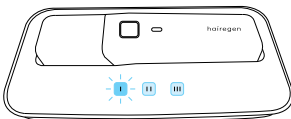STIMULATION PROFILES

| Profile | Most suitable for | Duration | Required Pressure | Sessions per day |
|---------|--------------------------|----------|-------------------|------------------|
| I | Beginners | 3 min. | Light | 1-2 |
| II | Advanced users | 5 min. | Medium | 1-2 |
| III | Long or dense hair users | 8 min. | Medium | 1-2 |

* For the most updated information please visit our website: www.dimedico.co.il/hairegen


USING HAIREGEN

- Remove the applicator from the base.
- Shortly press the applicator's button. Within a few seconds, following a self-test period, the applicator emits a startup sound and begins to vibrate.
- To modify the vibration intensity, press the applicator's button repeatedly, cycling through 3 vibration intensity levels.
- To turn off the applicator, press the button continuously for a couple seconds, or place the applicator in the powered base. After a few seconds, the applicator emits a shutdown tune and turns off.
- The applicator automatically shuts down in any of the following events:
 - Applicator is cradled in a powered base.
 - Stimulation profile limit is reached.
 - Low battery level.
 - A malfunction occurs.
- For effective treatment, roll the applicator in short movements back and forth over the scalp areas requiring treatment. Repeat this motion over the entirety of the treated area until the session is completed.
- In case of intense use some discomfort and tingling from the micro currents may be felt. This actually indicates the treatment is vigorous and effective. If the discomfort level is high, reduce the pressure of the disks on the scalp.
- When there is suboptimal contact between the scalp and the disks, the applicator will beep. The cause may be insufficient pressure, having too few disks in contact with the scalp or too much hair under the disks. Try slowing down the rolling, improving the hold of the applicator and increasing the pressure of the device on the scalp.
- If the warning does not stop, your scalp might be too dry for Hairegen to work properly. If this happens, use the humidifier – fill with clean water and dampen your scalp, for long hair lightly dampen the partings. Do not wet your hair.
- When the contact with the scalp is sufficient, the applicator emits red laser light (LLLT) from between the disks. The red laser light turns off when contact is poor*.
- Beginners are encouraged to use stimulation profile I, and proceed to stimulation profile II a couple of weeks later. If the intensity and duration are too demanding, remain on profile I for longer.
- Profile III is intended for users with long or dense hair. For such users, it is recommended to part the hair and roll the applicator back and forth across the parting line in short movements for 5-10 seconds, then move on to another area or parting line, repeating the action for the duration of the treatment.

The applicator will shut down upon profile program completion.

REPLACING THE CARTRIDGE


The disks of the cartridge come in contact with the scalp, and thus the cartridge is for personal use only. With time and usage, the performance degrades. The usage of the cartridge is therefore monitored and the user is alerted when it needs to be replaced.

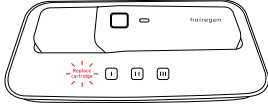When the cartridge usage exceeds 80% of its designed lifetime, a "Replace cartridge" alert will flash on the base. When the cartridge reaches its designed lifetime, the warning light "Replace cartridge" will remain on continuously on the base and the device will no longer turn on with that cartridge. When attempting to turn on the applicator with an empty cartridge, the light on the applicator will flash red a few times.

REMOVING THE CARTRIDGE

Hold the applicator with one hand and the protrusions at the cartridge's end with the other. Pull the cartridge out firmly. Push the new cartridge in until it clicks into place.

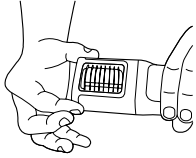MAINTENANCE

Keep all parts of the device clean and dry. The base and applicator should be cleaned with a dry cloth. Do not use solvents to clean any part of the Hairegen kit. The disks and the area between them can be cleaned with a brush. For sanitary purposes, roll the disks on a 70% alcohol pad, and dry them before use. Every 2 months, remove the cartridge and wipe the surfaces clean with a dry cloth.

DISPOSAL AND RECYCLING

The applicator contains a rechargeable battery. Please follow all relevant disposal regulations.

TECHNICAL SPECIFICATIONS

| Input | Output |
|---|---|
| <ul style="list-style-type: none">4.5 to 6 V DCPower consumption: up to 60mA when not charging, up to 500mA when charging.Interface: USB2 | <ul style="list-style-type: none">Voltage: up to 30VCurrent: up to 300µA, 0~100HzLaser: up to 5mW, 650 nm |

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